(DE) (BE)	TENS Pad Gebrauchsanweisung.....	4
(EN)	TENS Pad Instructions for use.....	14
(FR) (BE)	Coussinets TENS Mode d'emploi	22
(PL)	Elektrody TENS Pad Instrukcja użytkowania.....	32
(NL) (BE)	TENS Pad Gebruiksaanwijzing.....	41
(CZ)	TENS Pad Návod k použití	50
(SK)	TENS Pad Návod na použitie.....	58

Service-Hotline

(DE) Tel.: 0800 724 2355

(GB) Tel.: 0800 931 0319

(FR) Tél.: 0805 340 006


(NL) Tel.: 0800 543 0543

(BE) Tél.: 0800 706 11







(CZ) Tel.: 0800 555 013

(SK) Tel.: 0800 606 018

(PL) Tel.: 800 707 009

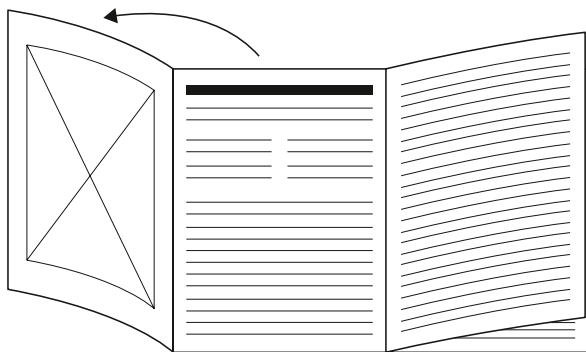- DE** Klappen Sie vor dem Lesen der Gebrauchsanweisung die Seite 3 aus.
- EN** Fold out page 3 before reading the instructions for use.
- FR** Dépliez la page 3 avant de lire le mode d'emploi.
- PL** Przed przeczytaniem instrukcji obsługi otworzyć stronę 3.
- NL** Vouw pagina 3 uit voordat u de gebruiksaanwijzing gaat lezen.
- CZ** Před přečtením návodu k použití rozložte stranu 3.
- SK** Pred čítaním návodu na použitie vyklopte stranu 3.

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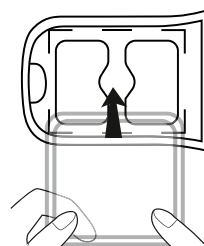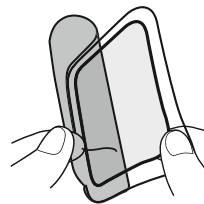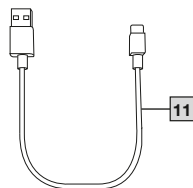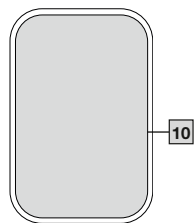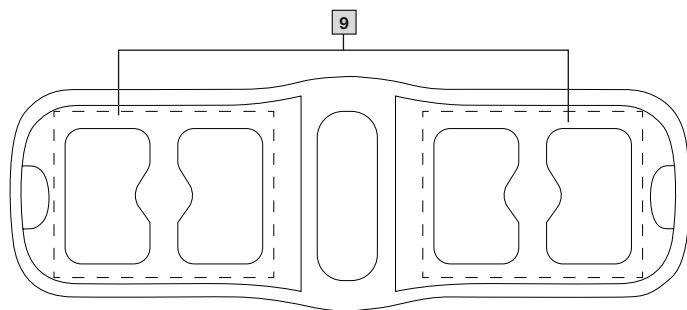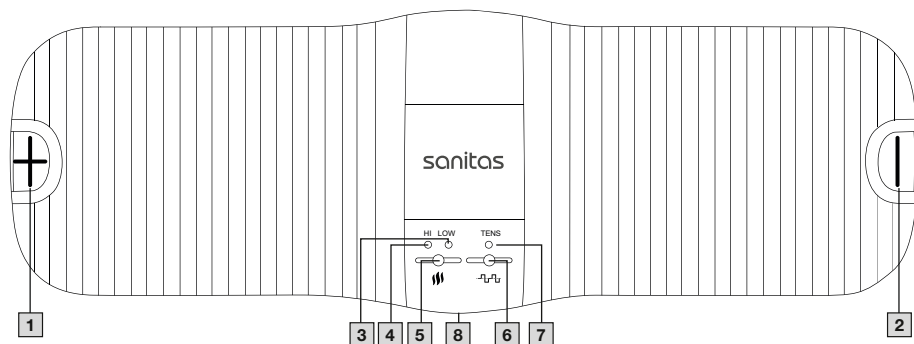Read these instructions for use carefully. Follow the warnings and safety notes. Keep these instructions for use for future reference. Make the instructions for use accessible to other users. If the device is passed on, provide the instructions for use to the next user as well.

CONTENTS

1. Included in delivery	14
2. Signs and symbols	14
3. Intended use	15
4. General warnings	16
5. Device description	18
6. Initial use	18
7. Usage	18
8. Cleaning and maintenance	19
9. Disposal	19
10. Accessories and replacement parts	20
11. Technical specifications	20
12. Guarantee/service	21

1. INCLUDED IN DELIVERY

Check that the exterior of the cardboard delivery packaging is intact and make sure that all contents are included in the delivery. Before use, ensure that there is no visible damage to the device or accessories and that all packaging material has been removed. If you have any doubts, do not use the device and contact your retailer or the specified Customer Service address.

- 1x SEM 53
- 1x USB charging cable
- 2x adhesive gel pads
- 1x instructions for use

2. SIGNS AND SYMBOLS

The following symbols are used on the device, in these instructions for use, on the packaging and on the type plate for the device:

⚠ WARNING

Indicates a potentially impending danger. If it is not avoided, there is a risk of death or serious injury.

⚠ CAUTION

Indicates a potentially impending danger. If it is not avoided, slight or minor injuries may result.

NOTICE

Indicates a potentially harmful situation. If it is not avoided, the device or something in its vicinity may be damaged.

Product information

Note on important information

Observe the instructions

Read the instructions before starting work and/or operating devices or machines

The device may not be used by people with medical implants (e.g. pacemakers), as this may affect their functionality.

Temperature range

Humidity range

Atmospheric pressure limitation

Use by

	Separate the product and packaging elements and dispose of them in accordance with local regulations.
	The device can emit effective output values above 10mA, averaged over every 5-second interval
	This product satisfies the requirements of the applicable European and national directives.
	Applied part, type BF
	Manufacturer
	Date of manufacture
	The electronic device must not be disposed of with household waste
	Do not dispose of batteries containing harmful substances with household waste
	Device protected against foreign objects ≥ 12.5 mm and against water dripping at an angle
	Serial number
	Marking to identify the packaging material. A = material abbreviation, B = material number: 1-7 = plastics, 20-22 = paper and cardboard
	Importer symbol
	Medical device
	Article number
	Type number
	Batch designation
	Authorised representative in the European Community

The device is Magnetic Resonance (MR) unsafe with a projectile Hazard)

3. INTENDED USE

The device is only intended for the purpose described in these instructions for use. The manufacturer is not liable for damage resulting from improper or careless use. We accept no liability for damage to the product, damage to property or injury to persons due to careless, improper or incorrect use or use of the product not in accordance with the purpose specified by the manufacturer. This device is intended for use in adults for temporary relief of back pain, joint pain, neuralgia, pain following musculoskeletal injury, pain due to circulatory disorders and chronic pain of various causes.


Important safety notes

Use of the device is not a substitute for medical consultation and treatment. Therefore, always consult your doctor first.

⚠ WARNING!

In order to avoid damage to health, do not use a stimulation current device in the following cases or mentioned body parts:

- If you have implanted electrical devices (e.g. a pacemaker)
- If you have metal implants
- People who use an insulin pump
- If you have a high temperature (e.g. $>39^{\circ}\text{C}$)
- If you have a known or acute cardiac arrhythmia or disorders of the heart's impulse and conduction system
- If you suffer from a seizure disorder (e.g. epilepsy)
- If you are pregnant
- If you have cancer
- After an operation, if strong muscle contractions could affect the healing process
- The device must never be used near to the heart. The stimulation electrodes must not be placed on any part of the front ribcage (where the ribs and breastbone are located), and especially not on the two large pectorals. This can increase the risk of ventricular fibrillation and induce cardiac arrest.
- On the skeletal skull structure, or around the mouth, throat or larynx
- In the neck/carotid artery area

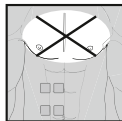- Do not use on the genital area
- On acutely or chronically diseased (injured or irritated) skin (e.g. inflamed skin – whether painful or not, reddened skin, rashes, e.g. allergies, burns, bruises, swellings, both open and healing wounds, and post-operative scars where the healing process could be affected)
- In humid environments (e.g. in the bathroom) or when bathing or showering
- Do not use after consuming alcohol
- If you are connected to a high-frequency surgical device
- Do not use if you have an acute or chronic disease of the gastrointestinal tract
- The stimulation must not be applied above or through the head, directly onto the eyes, on the mouth, to the front of the neck (especially not to the carotid artery), or with the electrodes or adhesive pads placed on the chest and upper back or across the heart
- The device does not be used in Magnetic Resonance (MR) environment because it is unsafe and with a projectile hazard.

Before using the device, consult your doctor if any of the following apply to you:

- If you suffer from an acute medical condition, in particular if you suspect or have been diagnosed with high blood pressure, a blood coagulation disorder, propensity to thrombo-embolic conditions or recurrent malignant growths
- If you have any skin conditions
- If you have unexplained chronic pain in any part of the body
- If you suffer from diabetes
- If you have any sensory impairment that reduces the feeling of pain (e.g. metabolic disorders)
- If you are receiving medical treatment
- In the event of complaints linked to the stimulation treatment
- If you suffer from persistently irritated skin due to long-term stimulation at the same electrode site

⚠ WARNING!

Only use this Digital TENS/EMS device:

- on adults.
- for the intended purpose and as specified in these instructions for use. Any form of improper use can be dangerous.
- for external use.

- with the original accessories supplied, which can be re-ordered. Failure to do so invalidates the guarantee.

4. GENERAL WARNINGS

- Always pull gently on adhesive gel pads to remove them from the skin to prevent injuries in the unusual case of highly sensitive skin.
- Hold the device away from sources of heat and do not use it in close proximity (approx. 1 m) to shortwave or microwave devices (e.g. mobile phones), as doing so can result in unpleasant current peaks.
- Do not expose the device to direct sunlight or high temperatures.
- Protect the device from dust, dirt and moisture.
- Never immerse the device in water or other liquids.
- The device is suitable for self-treatment.
- For hygiene reasons, the adhesive gel pads may only be used on one person.
- If the device does not work properly, or if you feel unwell or experience pain, stop using it immediately.
- Switch off the device first before removing or relocating the device to prevent unintentional stimulation.
- Do not modify electrodes or adhesive gel pads (e.g. by cutting them). Doing so increases the current density, which is potentially hazardous (max. recommended output value for the electrodes is 9 mA/cm², an effective current density beyond 2 mA/cm² requires increased caution).
- Do not use whilst asleep, driving a vehicle or operating machinery.
- Do not use whilst undertaking any activity where an unexpected reaction (e.g. strong muscle contractions even at low intensity) could be dangerous.
- Ensure that no metallic objects (e.g. belt buckles or necklaces) come into contact with the electrodes during stimulation. If you are wearing jewellery or have piercings in the area to be treated (e.g. a navel piercing), these must be removed before using the device; failure to do this could result in spot burns.
- Keep the device away from children to prevent potential hazards.
- Make sure not to confuse the electrode cables and their contacts with your headphones or other devices, and do not connect the electrodes to other devices.
- Do not use this device whilst using other devices that transmit electrical pulses to your body.
- Do not use in the vicinity of highly flammable substances, gases or explosives.

- During the initial few minutes, use the device while sitting or lying down to minimise the risk of accidental injuries as a consequence of isolated cases of vagal responses (feeling of faintness). If you feel faint, switch off the device immediately, lie down and support the legs in an elevated position (approx. 5–10 minutes). Treatment of the skin with moisturising lotions or ointments beforehand is not recommended as this considerably increases the electrode wear and may cause unpleasant current peaks.
- This device is not intended for use by children or people with restricted physical, sensory (e.g. reduced sensitivity to pain) or mental skills, or people lacking experience and/or knowledge, unless they are supervised by a person who is responsible for their safety or are instructed by such a person in how to use the device.
- If the adhesive capability of the adhesive gel pads decreases, please replace them immediately. Do not use the device until you have fitted the new adhesive gel pads. Otherwise the uneven adhesion of the adhesive gel pads could lead to skin injuries.
- Replace the adhesive gel pads with new pads at the latest after having used them 20 times.

Damage

- If you notice any signs of wear or damage, do not use the device and contact your retailer or the specified Customer Services address.
- Switch the device off immediately if it is faulty or not working properly. Never attempt to open or repair the device yourself. Repairs may only be carried out by the Service Centre or authorised specialists. Failure to comply with this instruction will void the guarantee.
- The manufacturer is not liable for damage resulting from improper or incorrect use.

Notes on handling batteries

WARNING

- **Risk of explosion! Risk of fire!** Failure to comply with the following points can result in personal injury or cause overheating, leakage, venting, breakage, explosion, or fire on the battery.
- This device contains non-rechargeable batteries which must not be charged.
- Do not throw batteries into a fire.
- Never charge, forcibly discharge, heat, disassemble, open, crush, deform, encapsulate, or modify batteries.
- Never short-circuit batteries or battery compartment contacts.

- Protect the batteries from direct sunlight, rain, heat, and water.
- Exposure of batteries to an environment with extremely high temperatures or an extremely low air pressure may result in explosion or leakage of flammable liquids and gases.
- Dispose of defective and discharged batteries immediately and properly (see chapter on disposal).
- Do not use modified or damaged batteries.
- Always select the correct battery type.
- Always insert the batteries correctly, taking into account the polarity (+ / -).
- Never mix batteries of different manufacturers, capacities (new and used), size, or type within a device.
- If a battery has leaked, put on protective gloves and clean the battery compartment with a dry cloth.
- If fluid from a battery comes into contact with your skin or eyes, wash the affected areas with water and seek medical assistance.
- **Choking hazard!** Keep batteries out of the reach of children. Seek medical attention immediately if swallowed. Swallowing them may cause burns, severe internal injuries, and death.
- Never allow children to replace batteries without adult supervision.

CAUTION

- Store batteries in a well-ventilated, dry, and cool place in a non-conductive container in which the batteries cannot be short-circuited to each other or by other metal objects.
- Keep batteries clean and dry.
- Keep batteries away from water.
- If the device is not going to be used for a long period of time, remove the batteries from the battery compartment.
- In order to achieve a battery service life as long as possible, fully charge the battery at least every 3–6 months

NOTICE

- Do not use rechargeable batteries.

This device can be used to combat pain (TENS = transcutaneous electrical nerve stimulation).

TENS stands for transcutaneous (through the skin) electrical nerve stimulation. Chronic pain can be treated with various methods, e.g. with medications, by

surgery, but also with other methods without side effects, such as acupuncture or TENS.

5. DEVICE DESCRIPTION

- 1 Intensity button+, button lock
- 2 Intensity button -
- 3 Heat setting indicator LED, low
- 4 Heat setting indicator LED, high
- 5 Heat function button HEAT
- 6 Function/pause button TENS
- 7 TENS indicator LED
- 8 USB-C charging socket
- 9 Electrodes
- 10 Adhesive gel pads (2x)
- 11 USB cable

6. INITIAL USE

Charging the battery

 Allow the device to charge for at least 2 hours before first use.

- Connect the supplied USB cable to a USB power source and the USB-C socket on the device.
- Plug the USB charger (SV/1 A) into a suitable socket.
- The TENS indicator LED lights up green during the charging process. When the battery is fully charged, the TENS indicator LED goes out.

You cannot use this device while it is charging. Only use this product with a suitable USB cable and a suitable USB power source that meets the specially applicable standards.

Applying the gel pads

- Remove one of the protective films carefully .
- Apply the gel pads to the electrode area, ensuring that the edge of the gel pad does not protrude beyond the electrode profile. Apply only one gel pad on each side so that it covers the two electrodes on one side .
- Remove the second protective film slowly and carefully to ensure that the gel pads are not damaged.

NOTICE

Only use this product with the original gel pads supplied, which can be reordered. Failure to follow this rule shall void the guarantee. The manufacturer is not liable for damage resulting from improper or careless use.

7. USAGE

- There are 3 modes to choose from (TENS mode, heat mode and combined TENS/heat mode).
- This device has two heat settings and 15 intensity settings for TENS use.
- We recommend that you start with a low intensity level and gradually increase it after getting used to the device.
- This device has an automatic switch-off function. It switches off automatically after 30 minutes of combined TENS and heating use or after 20 minutes of heating only or TENS only use.

Step 1: Please adhere the device to the painful area using a gel pad. This device has a function that detects whether the device is correctly positioned on the body. If there is not sufficient contact with the body, the device beeps and switches itself off.

Step 2: Select one of the following modes:

A) TENS mode:

- Press the TENS button  or the "+" button  for 1 second to switch on the device and start the TENS function. The TENS display  lights up green.
- Briefly press the "+" button  to increase the intensity.
- Briefly press the "-" button  on the device to reduce the intensity.
- Button lock: Press and hold the "+" button  for 2 seconds to activate or deactivate the button lock. An acoustic signal sounds when the button

lock is activated and deactivated. If the button lock is activated, the TENS intensity or heat setting cannot be adjusted, but the device can be switched off by pressing and holding the “-” button [2].

- If necessary, briefly press the TENS button [6] to pause or resume TENS mode.
- If necessary, briefly press the HEAT button [5] to switch to the combined TENS and heat mode (TENS and heat indicator LEDs [6] [3] or [4] light up).

B) Heat mode

- Press the HEAT button [5] for 1 second to start the low heat setting. The corresponding indicator LED [3] lights up green.
- Press the HEAT button [5] again to switch to the high heat setting. The corresponding indicator LED [4] lights up green.
- **Button lock:** Press and hold the “+” button [1] for 2 seconds to activate or deactivate the button lock. An acoustic signal sounds when the button lock is activated and deactivated. If the button lock is activated, the TENS intensity or heat setting cannot be adjusted, but the device can be switched off by pressing and holding the “-” button [2].
- If necessary, briefly press the TENS button [6] to switch to the combined TENS and heat mode (TENS and heat indicator LEDs [6] [3] or [4] light up).

C) Combined TENS and heat mode

- Press the TENS [6] and HEAT [5] buttons simultaneously to activate the combined mode. Alternatively, press the TENS button first and then the HEAT button [5]. The corresponding indicator LED lights up green.
- Briefly press the “+” button [1] to increase the intensity. Briefly press the “-” button [2] on the device to reduce the intensity.
- **Button lock:** Press and hold the “+” button [1] for 2 seconds to activate or deactivate the button lock. An acoustic signal sounds when the button lock is activated and deactivated.
- If the button lock is activated, the TENS intensity or heat setting cannot be adjusted, but the device can be switched off by pressing and holding the “-” button [2]. If necessary, briefly press the TENS button [6] to pause or resume TENS mode.
- If necessary, briefly press the HEAT button [5] again to switch from low to high heat setting. If the heat function is activated, the heat function is switched off by pressing and holding the HEAT button [5].

Step 3: As described above, the device switches off automatically after either 20 or 30 minutes of operation. The device can be switched off manually by pressing and holding the “-” button [2] or the TENS button [6].

This device has a function that detects whether the device is correctly positioned on the body. If there is not sufficient contact with the body, the device beeps and switches itself off. If the device is attached to your body but still beeps, clean the gel pads and try starting the device again. If the continuous beeping persists, replace the gel pads.

Usage instructions

- Switch off the device before removing it.
- Wait at least 30 minutes between applications.
- Only use the device as often as you find comfortable. If the treatment with this device feels uncomfortable, stop using it.
- If the gel pads lose their adhesion, replace them immediately. Do not use the device again until the gel pads have been replaced, as uneven adhesion of the gel pads may lead to skin lesions.
- In the event of unpleasant current peaks, switch off the device and slightly change its position.

8. CLEANING AND MAINTENANCE

Never use solvents or abrasives, hard brushes, metal or sharp objects for cleaning. Solvents are harmful to health and weaken plastic parts; abrasive cleaners and tools scratch the surface.

- Clean the outer surfaces of the product with a slightly damp cloth and then dry thoroughly.
- Clean the adhesive gel pads with a damp, lint-free cloth or under lukewarm, running water and pat dry with a lint-free cloth.
- Store the product in a dry, cool place away from moisture and direct sunlight. Store the product out of the reach of children.

9. DISPOSAL

For environmental reasons, do not dispose of the device in household waste at the end of its service life. Dispose of the device at a suitable local collection or recycling point in your country. Observe the local regulations for material disposal. Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal. You can obtain the location of collection points for old devices from the local or municipal authorities, local waste disposal companies or your retailer, for example.

Empty, completely flat batteries must be disposed of through specially designated collection boxes, recycling points, or at electronics retailers. You are legally required to dispose of the batteries appropriately.

The codes below are printed on batteries containing harmful substances:

Pb = battery contains lead

Cd = battery contains cadmium

Hg = battery contains mercury

10. ACCESSORIES AND REPLACEMENT PARTS

Accessories are available to purchase from the webshop at:

www.shop.sanitas-online.de.

Designation	Item number and/or order number
SEM 53 Gel Pads	110.344

11. TECHNICAL SPECIFICATIONS

Type	FT-810
Model	SEM 53
Output waveform	Biphasic rectangular pulses
Pulse duration	200 - 400 µs
Pulse frequency	6 - 45 Hz
Output voltage	Max. 80 Vp-p (at 1 kΩ)
Output current	80 mA _p -p (at 1 kΩ)
Heat mode	HI: up to 43°C LOW: up to 41°C
Nominal voltage	3.7 V — — —
Nominal input voltage/current	5 V — — —, 1 A
Battery	Li-Ion, 3.7 V, 1000 mAh
Treatment time, automatic switch-off	TENS mode: 20 minutes, TENS and heat mode: 30 minutes, Heat mode: 20 minutes
Intensity levels	0 to 15

Operating conditions:	5-40°C at a relative humidity of 15-90%, atmospheric pressure: 700 hPa to 1060 hPa
Storage conditions	0-40°C at a relative humidity of 0-90%
Dimensions	100 mm x 300 mm x 15 mm
Weight	89 g
Gel pad usage cycles	20
Expected service life of the device	Information on the service life of the product can be found on the website

This device complies with the European standard EN 60601-1-2 (Group 1, Class B, in compliance with CISPR 11, IEC 61000-4-2, IEC 61000-4-3, IEC 61000-4-8) and is subject to special precautionary measures with regard to electromagnetic compatibility. Please note that portable and mobile HF communication systems may interfere with this device. For more details, please contact Customer Services at the address indicated.

This device meets the requirements of Directive 93/42/EEC for medical devices, as well as those of the Medizinproduktegesetz (German Medical Devices Act), as applicable. For this device, no additional operator-performed functional test and user instruction are required, and no safety checks are required under the German Medical Device Operator Ordinance (MPBetreibV) provisions on user instruction/functional checks and safety-related inspections, as applicable.

Notes on electromagnetic compatibility

- The device is suitable for use in all environments listed in these instructions for use, including domestic environments.
- The use of the device may be limited in the presence of electromagnetic disturbances. This could result in issues such as error messages or the failure of the display/device.
- Avoid using this device directly next to other devices or stacked on top of other devices, as this could lead to faulty operation. If, however, it is necessary to use the device in the manner outlined above, this device as well as the other devices must be monitored to ensure they are working properly.
- The use of accessories other than those specified or provided by the manufacturer of this device can lead to an increase in electromagnetic emissions or a decrease in the device's resistance to electromagnetic interference; this can result in faulty operation.

- Keep portable RF communication devices (including peripheral equipment, such as antenna cables or external antennas) at least 30 cm away from all device parts, including all cables included in delivery.
- Failure to comply with the above could impair the performance of the device.

12. GUARANTEE/SERVICE

Hans Dinslage GmbH, Riedlinger Straße 28, 88524 Uttenweiler (hereinafter referred to as "HaDi"), provides a guarantee for this product, subject to the requirements below and to the extent described as follows.

Claiming the guarantee shall not affect the Purchaser's statutory guarantee rights arising from the purchase agreement with the Seller in the event of defects. The Purchaser may exercise these statutory guarantee rights free of charge. The guarantee shall apply without prejudice to any mandatory statutory provisions on liability.

HaDi guarantees the perfect functionality and completeness of this product.

The worldwide guarantee period is 3 years, commencing from the purchase of the new, unused product from the seller.

The guarantee only applies to products purchased by the buyer as a consumer and used exclusively for personal purposes in the context of domestic use. German law shall apply.

During the guarantee period, should this product prove to be incomplete or defective in functionality in accordance with the following provisions, HaDi shall carry out a repair or a replacement delivery free of charge, in accordance with these guarantee conditions.

Should the buyer wish to file a guarantee claim, they shall first contact the HaDi Customer Service team:

Service hotline (free):

☎ Tel.: 0800 931 0319

E-mail:

service-uk@sanitas-online.de

The buyer will then receive further information about the processing of the guarantee claim, e.g. where they can send the product and what documentation is required.

If we request the buyer to send the defective product, the product must be sent to the following address:

NU Service GmbH, Lessingstraße 10 b
89231 Neu-Ulm, Germany

A guarantee claim shall only be considered if the buyer can provide HaDi, or an authorised HaDi partner, with

- a copy of the invoice/purchase receipt, and
- the original product.

The following are explicitly excluded from this guarantee:

- deterioration due to normal use or consumption of the product;
- accessories supplied with this product which are worn out or used up through proper use (e.g. batteries, rechargeable batteries, cuffs, seals, electrodes, light sources, attachments and nebuliser accessories);
- products that are used, cleaned, stored or maintained improperly and/or contrary to the provisions of the instructions for use, as well as products that have been opened, repaired or modified by the buyer or by a service centre not authorised by HaDi;
- damage that arises during transport between manufacturer and customer, or between service centre and customer;
- products purchased as seconds or as used goods;
- consequential damage arising from a fault in this product (however, in this case, claims may exist arising from product liability or other compulsory statutory liability provisions).

Repairs or an exchange in full do not extend the guarantee period under any circumstances.

Notification of incidents

For users/patients in the European Union and identical regulation systems, the following applies: If a major incident occurs during or through use of the product, notify the manufacturer and/or their representative of this as well as the respective national authority of the member state in which the user/patient is located.

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HIVOX BIOTEK B.V.