sanitas

SBR 05

(DE) (BE) Insektenstichheiler Gebrauchsanweisung	4
(GB) Insect bite healer Instructions for use	37
(FR) (BE) Appareil soulageant les piqûres d'insectes Mode d'emploi.....	65
(NL) (BE) Insectenbeetheler Gebruiksaanwijzing	99

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Read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

TABLE OF CONTENTS

1. Included in delivery	38	8. Cleaning and maintenance	55
2. Signs and symbols	38	9. Replacement parts	56
3. Proper use	42	10. What if there are problems?	57
4. Warnings and safety notes	44	11. Disposal	59
5. Device description	51	12. Technical specifications	59
6. Initial use	51	13. Warranty/service	61
7. Usage	52		

1. INCLUDED IN DELIVERY

Check that the exterior of the cardboard delivery packaging is intact and make sure that all contents are present. Before use, ensure that there is no visible damage to the device and karabiner and that all packaging material has been removed. In case of any doubt, do not use the device, and contact your retailer or the specified Customer Service address.

1 x SBR 05

1 x Carabiner

2 x batteries, type AAA (LR03 1.5 V)

1 x instructions for use

2. SIGNS AND SYMBOLS

The following symbols are used on the device, in these instructions for use, on the packaging and on the type plate for the device:

WARNING

Indicates a potentially impending danger. If it is not avoided, there is a risk of death or serious injury.

CAUTION

Indicates a potentially impending danger. If it is not avoided, slight or minor injuries may result.

Product information

Note on important information

Observe the instructions

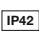
Read the instructions before starting work and/or operating devices or machines

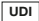
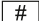
Type BF applied part

On

○	Off
	The electronic device must not be disposed of with household waste
	Do not dispose of batteries containing harmful substances with household waste
	Manufacturer
	Serial number
	Device protected against solid foreign objects ≥ 1 mm in diameter and larger, and against drops of water when the housing is angled up to 15°
	Temperature limit
	Humidity limitation

	<p>CE labelling This product satisfies the requirements of the applicable European and national directives.</p>
	<p>Medical device</p>
	<p>Extra-small design which is particularly suitable for use on the move</p>
	<p>Separate the product and packaging elements and dispose of them in accordance with local regulations.</p>
	<p>Marking to identify the packaging material. A = Material code, B = Material number: 1-7 = Plastics, 20-22 = Paper and cardboard</p>
	<p>Unique Device Identifier (UDI) for unique product identification</p>
	<p>Item number</p>
	<p>Type number</p>
	<p>Date of manufacture</p>

LOT	Batch code
------------	------------

3. PROPER USE

Intended use

The insect bite healer is intended for relieving itching and swelling caused by insect bites or stings by applying heat (local hyperthermia).

Target group

The insect bite healer is intended for use in the home healthcare environment, not in professional healthcare facilities. The insect bite healer can be used on anyone over 3 years of age under supervision; can be used for self-treatment by anyone over 12 years of age.

Clinical benefits

Treatment and relief of itching and swelling caused by insect bites or stings.

Indication

Treatment with the insect bite healer is recommended for:

- itching

- swelling
caused by insect bites or stings.

CONTRAINDICATIONS

WARNING

- Do not use the device if you suffer from diabetes. Sensitivity to pain could, in certain circumstances, be reduced and diabetics may burn themselves.
- Do not use the device on people with disabilities, children under 3 years old or people with reduced sensitivity to heat. This includes people with skin alterations due to illness or scarred tissue in the application area, after taking pain relief medication, drinking alcohol or taking drugs.
- Do not use the device on acutely or chronically diseased (injured or irritated) skin (e.g. inflamed skin, whether painful or not, reddened skin, rashes, (e.g. allergies), burns, frostbite, bruises, swellings (without signs of an insect bite/sting), open/bleeding and healing wounds, and surgical scars where the healing process could be affected).
- Do not use the device if you are suffering from a fever.

- Do not use the device on sensitive areas of skin.
- Do not use the device if you have any sensory impairment that reduces the feeling of pain (e.g. metabolic disorders).
- Do not use the device if you are using lotions, creams and gels locally at the same time.
- Do not use the device if you suffer from persistently irritated skin due to long-term heat application on the same area of skin.

4. WARNINGS AND SAFETY NOTES

GENERAL WARNINGS

⚠ WARNING

- Any person using it on themselves must be over 12 years old. The product should only be used on children under 12 years old if supervised or assisted by an adult.
- Before using the device, consult your doctor if you have any skin conditions, serious illness (in particular if you have a propensity to thromboembolic con-

ditions or recurrent malignant growths), unexplained chronic pain in any part of the body.

- Seek medical advice immediately if the symptoms persist despite using the insect bite healer and stop using it.
- Seek medical advice immediately if the symptoms worsen whilst or after using the insect bite healer and stop using it.
- Seek medical advice immediately if it is a tick bite. The bite may result in transmission of pathogens such as tick-borne encephalitis (TBE) or borreliosis.
- Seek medical advice immediately if the first signs of increased itching, skin redness, swelling, dizziness, nausea, shortness of breath, formation of wheals or a drop or increase in blood pressure occur.
- Do not use near to the eyes, on the eyelids, mouth or on mucous membranes.
- Do not allow children to play with the device.
- Store the device out of the reach of children and pets.
- Keep packaging material away from children (risk of suffocation).
- Please note the general hygiene measures when using the insect bite healer.
- The device is not a substitute for medical consultation and treatment. Always consult your doctor if you are experiencing any pain or are suffering from an illness.

- If the device does not work properly, or if you feel unwell or experience pain, stop using it immediately.

GENERAL PRECAUTIONS

⚠ CAUTION

- Do not touch the device with wet hands. Do not allow water to spray onto the device. Only operate the device if it is completely dry.
- Do not use any detergents or solvents to clean and maintain the device. If liquid penetrates the device this could cause damage to the electronics or other device parts and lead to a malfunction.
- The device should only be used by one person; use by multiple people is not recommended.
- Protect the device against direct sunlight and store in a dry place.
- Protect the device from impacts, moisture, dust, chemicals, marked temperature fluctuations and nearby sources of heat (ovens, heaters).
- Ensure that no water enters the device. If this should happen, remove the batteries and only use the device again when it is completely dry.

- If you store the device, no heavy objects should be placed on top of it.
- Before use, ensure that there is no visible damage to the device and karabiner. In case of any doubt, do not use the device, and contact your retailer or the specified Customer Service address.
- Store the device in a location protected against climatic influences. The device must be stored in the environmental conditions specified.
- Do not leave the device unattended while in use. If the device malfunctions, switch it off and stop using it.

MEASURES FOR HANDLING BATTERIES

WARNING

- **Risk of explosion! Risk of fire!** Failure to comply with the following points can result in personal injury or cause overheating, leakage, venting, breakage, explosion, or fire on the battery.
- This device contains non-rechargeable batteries which must not be charged.
- Do not throw batteries into a fire.

- Never charge, forcibly discharge, heat, disassemble, open, crush, deform, encapsulate, or modify batteries.
- Never short-circuit batteries or battery compartment contacts.
- Protect the batteries from direct sunlight, rain, heat, and water.
- Exposure of batteries to an environment with extremely high temperatures or an extremely low air pressure may result in explosion or leakage of flammable liquids and gases.
- Dispose of defective and discharged batteries immediately and properly (see chapter on disposal).
- Do not use modified or damaged batteries.
- Always select the correct battery type.
- Always insert the batteries correctly, taking into account the polarity (+ / -).
- Never mix batteries of different manufacturers, capacities (new and used), size, or type within a device.
- If a battery has leaked, put on protective gloves and clean the battery compartment with a dry cloth.
- If fluid from a battery comes into contact with your skin or eyes, wash the affected areas with water and seek medical assistance.

- **Choking hazard!** Keep batteries out of the reach of children. Seek medical attention immediately if swallowed. Swallowing them may cause burns, severe internal injuries, and death.
- Never allow children to replace batteries without adult supervision.

⚠ CAUTION

- Store batteries in a well-ventilated, dry, and cool place in a non-conductive container in which the batteries cannot be short-circuited to each other or by other metal objects.
- Keep batteries clean and dry.
- Keep batteries away from water.
- If the device is not going to be used for a long period of time, remove the batteries from the battery compartment.

NOTICE

- Do not use rechargeable batteries.

NOTE ON ELECTROMAGNETIC COMPATIBILITY

WARNING

- The device is suitable for use in all environments listed in these instructions for use, including domestic environments.
- The use of the device may be limited in the presence of electromagnetic disturbances. This could result in issues such as error messages or the failure of the device.
- Avoid using this device directly next to other devices or stacked on top of other devices, as this could lead to faulty operation. If, however, it is necessary to use the device in the manner stated above, this device as well as the other devices must be monitored to ensure they are working properly.
- The use of accessories other than those specified or provided by the manufacturer of this device can lead to an increase in electromagnetic emissions or a decrease in the device's electromagnetic immunity; this can result in faulty operation.

- Keep portable RF communication devices (including peripheral equipment, such as antenna cables or external antennas) at least 30 cm away from all device parts, including all cables included in delivery.
- Failure to comply with the above can impair the performance of the device.

5. DEVICE DESCRIPTION

The associated drawings are shown on page 3.

1 Activation button

4 Heating plate

2 ON/OFF switch

5 Battery compartment lid

3 Device status LED

6 Carabiner

6. INITIAL USE

Inserting the batteries

1. Push open the battery compartment lid **5**.

2. Insert the batteries included in delivery into the battery compartment as shown inside the battery compartment. When inserting the batteries, ensure that the correct polarity (+/-) is observed.
3. Slide the battery compartment lid **5** closed until you feel and hear it click into place.

7. USAGE

⚠ WARNING

Before use, make sure that there is no insect sting remaining in the area that was stung. Remove the sting carefully before use. To achieve the best results, use the device immediately after the insect bite or sting. If you wait too long to treat the area that has been stung or bitten and itching and swelling have already started, the insect bite healer can only have a limited effect. Nevertheless, treat the area with the insect bite healer. Usually, itching can be relieved and healing can be sped up.

CAUTION

Depending on the sensitivity of the person using the device and of the area of skin, the treatment temperature may feel unpleasant and cause slight skin redness and, in rare cases, irritated skin. If the heat function feels too hot, stop treatment immediately!

1. To switch on the device, slide the ON/OFF control **2** to the “I” position in the direction of the heating plate **4**.
→ Device status LED **3** starts to light up blue. Device is now ready for use.
2. Place the heating plate **4** of the device on the area that has been stung or bitten.

The 3-second treatment is recommended for first use or use on sensitive skin.

- To use the 3-second treatment, **briefly** press the activation button **1** once.
→ The device status LED **3** lights up blue and the treatment starts.

If you press and hold the activation button **1** for more than 2 seconds, the emergency switch-off will be activated. To use the device again after the automatic switch-off, switch the device back to “OFF” using the ON/OFF control **2**. This is indicated by the “O” symbol. The blue device status LED **3** turns off.

3. An acoustic signal sounds at the end of the treatment.

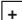
→ The device status LED **3** lights up green again and the heating plate **4** stops providing heat.

4. To prolong the battery life, switch off the device each time after using it. To switch off the device, slide the ON/OFF control **2** to the “O” position.

If you do not actively switch off the device yourself, the automatic switch-off function will switch it off automatically after 2 minutes. To use the device again after the automatic switch-off, switch the device back to “OFF” using the ON/OFF control **2**. This is indicated by the “O” symbol. The blue device status LED **3** turns off.

Battery warning

Before the battery level becomes too low for other applications, the device status LED  starts to flash 3x and an acoustic signal sounds 3x. Replace both batteries soon, as described in section 6.

8. CLEANING AND MAINTENANCE

Observe the following hygiene instructions to avoid health risks:

- The device is intended for multiple uses.
- It is recommended that the device is cleaned before first use and each additional use.
- The device must be switched off, allowed to cool down and the batteries removed each time before cleaning.
- Use only soft, dry cloths for cleaning the housing.
- Never use any substances for cleaning that could potentially be toxic if they came into contact with the skin or mucous membranes, or if they were swallowed or inhaled. Do not use any abrasive cleaning products and never submerge the device in water.
- Do not clean the device and karabiner in the dishwasher!

- Avoid mechanical cleaning with brushes or similar, as this could result in irreparable damage.
- When using the device on multiple areas that have been stung and bitten, it is recommended that you clean the heating plate each time. This reduces the possibility of germs and bacteria being transferred from one area of skin to the next. To disinfect the heating plate, use a cotton bud moistened with disinfectant.

9. REPLACEMENT PARTS

The replacement part is available from the corresponding listed service address under the stated material number.

Designation	Item or order number
Carabiner	164.395

10. WHAT IF THERE ARE PROBLEMS?

Problem	Possible cause	Solution
The device does not switch on.	The ON/OFF switch 2 is in the OFF position.	Slide the ON/OFF switch 2 into the ON position.
	The batteries are flat.	Replace both batteries as described in these instructions for use.
The device no longer gets hot enough.	The battery level is too low.	Replace both batteries as described in these instructions for use.
Device outputs 3 acoustic signals/the device status LED 3 flashes – even after switching the device OFF and ON again.	Battery level is too low, meaning that the device may no longer be able to be operated.	Replace both batteries as described in these instructions for use.

Problem	Possible cause	Solution
<p>Application cannot be started, even though the ON/OFF switch 2 is set to ON.</p>	<p>Device has switched off automatically. Either the activation button was pressed for more than 2 seconds or it was not used for more than 2 minutes and the automatic switch-off was activated.</p>	<p>First switch the device to OFF using the ON/OFF switch 2 and then switch it on again as usual.</p>
<p>The device does not switch on, even though the batteries have been replaced and the ON/OFF switch 2 has been set to ON.</p>	<p>The device may be faulty.</p>	<p>Please contact Customer Services.</p>

11. DISPOSAL

For environmental reasons, do not dispose of the device in household waste at the end of its service life. Dispose of the device at a suitable local collection or recycling point in your country. Dispose of the device in accordance with EC Directive Waste Electrical and Electronic Equipment (WEEE). If you have any questions, please contact the local authorities responsible for waste disposal.

Batteries must not be disposed of in the household waste. They may contain poisonous heavy metals and are subject to special refuse treatment. The codes below are printed on batteries containing harmful substances:

Pb = battery contains lead, Cd = battery contains cadmium

Hg = battery contains mercury

12. TECHNICAL SPECIFICATIONS

Type	SBR 55
Model	SBR 05

Dimensions (L x W x H)	86 x 35 x 25 mm
Weight (with batteries)	44.5 g
Power supply	2 x 1.5 V AAA LR03 batteries
Max. temperature	50 °C ±2 °C (122 °F ±3.6 °F)
Functions	1 time settings (3s)
Expected operating duration	Approx. 300 applications (per batteries)
Permissible storage and transport temperature and humidity	-10 °C to 50 °C, relative humidity 30% to 85%
Permissible operating temperature and humidity	10 °C to 40 °C, relative humidity 30% to 85%
Automatic device switch-off	2 min after last application

The serial number is located on the device or in the battery compartment.

Subject to technical changes.

The device complies the respective national regulations and the European standard EN 60601-1-2 (Group 1, Class B, in accordance with CISPR 11, CISPR 14-1, CISPR-32, IEC 61000-3-2, IEC 61000-3-3, IEC 61000-4-2, IEC 61000-4-3,

IEC 61000-4-4, IEC 61000-4-5, IEC 61000-4-6, IEC 61000-4-8, IEC 61000-4-11) and is subject to particular precautions with regard to electromagnetic compatibility.

13. WARRANTY/SERVICE

Hans Dinslage GmbH, Riedlinger Straße 28, 88524 Uttenweiler (hereinafter referred to as “HaDi”), provides a guarantee for this product, subject to the requirements below and to the extent described as follows.

Claiming the guarantee shall not affect the Purchaser's statutory guarantee rights arising from the purchase agreement with the Seller in the event of defects. The Purchaser may exercise these statutory guarantee rights free of charge. The guarantee shall apply without prejudice to any mandatory statutory provisions on liability.

HaDi guarantees the perfect functionality and completeness of this product.

The worldwide guarantee period is 3 years, commencing from the purchase of the new, unused product from the seller.

The guarantee only applies to products purchased by the buyer as a consumer and used exclusively for personal purposes in the context of domestic use. German law shall apply.

During the guarantee period, should this product prove to be incomplete or defective in functionality in accordance with the following provisions, HaDi shall carry out a repair or a replacement delivery free of charge, in accordance with these guarantee conditions.

Should the buyer wish to file a guarantee claim, they shall first contact the HaDi Customer Service team:

Service hotline (free):

 Tel.: 0800 931 0319

E-mail:

service-uk@sanitas-online.de

The buyer will then receive further information about the processing of the guarantee claim, e.g. where they can send the product and what documentation is required.

If we request the buyer to send the defective product, the product must be sent to the following address:

NU Service GmbH
Lessingstraße 10 b
89231 Neu-Ulm
Germany

A guarantee claim shall only be considered if the buyer can provide HaDi, or an authorised HaDi partner, with

- a copy of the invoice/purchase receipt, and
- the original product.

The following are explicitly excluded from this guarantee:

- deterioration due to normal use or consumption of the product;
- accessories supplied with this product which are worn out or used up through proper use (e.g. batteries, rechargeable batteries, cuffs, seals, electrodes, light sources, attachments and nebuliser accessories);
- products that are used, cleaned, stored or maintained improperly and/or contrary to the provisions of the instructions for use, as well as products that have been opened, repaired or modified by the buyer or by a service centre not authorised by HaDi;

- damage that arises during transport between manufacturer and customer, or between service centre and customer;
- products purchased as seconds or as used goods;
- consequential damage arising from a fault in this product (however, in this case, claims may exist arising from product liability or other compulsory statutory liability provisions).

Repairs or an exchange in full do not extend the guarantee period under any circumstances.

Notification of incidents

For users/patients in the European Union and identical regulation systems, the following applies: If during or through use of the product a major incident occurs, notify the manufacturer and/or their representative of this as well as the respective national authority of the member state in which the user/patient is located.

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